camry


CR 8414


(GB) user manual
(F) mode d'emploi
(P) manual de serviço
(LV) lietošanas instrukcija
(RO) Instrucțiunea de deservire
(HU) felhasználói kézikönyv
(MK) Упатство за корисникот
(RU) Инструкция по эксплуатации

(D) bedienungsanweisung
(E) manual de uso
(LT) naudojimo instrukcija
(EST) kasutusjuhend
(BIH) upute za uporabu
(CZ) návod k obsluze
(PL) instrukcja obsługi


A.

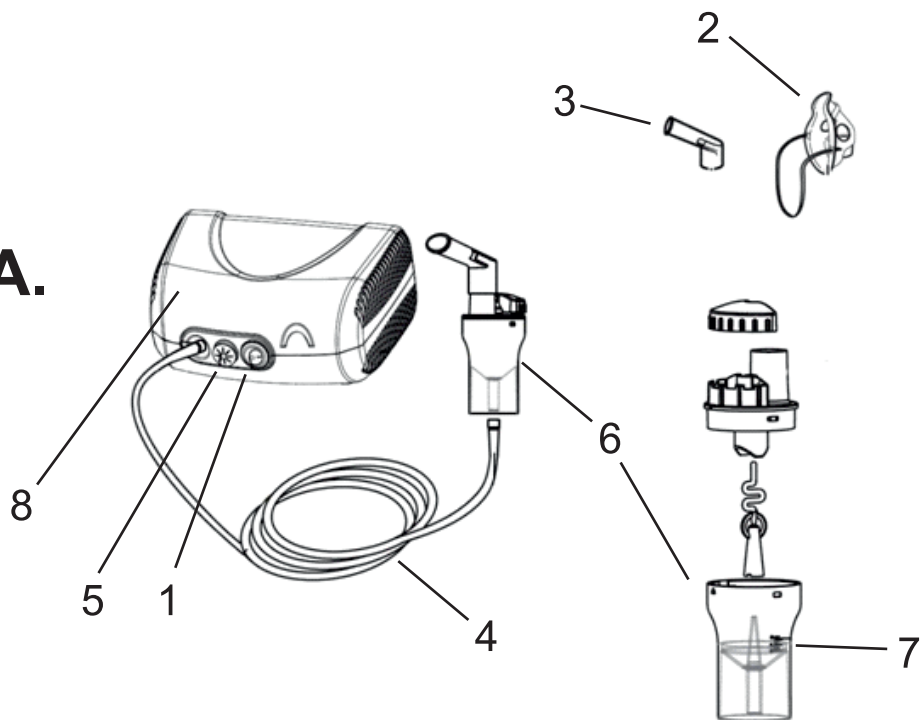
B.

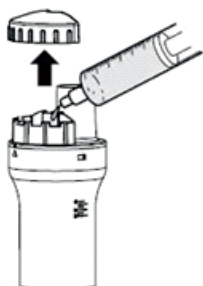C.

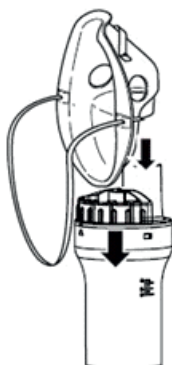D.

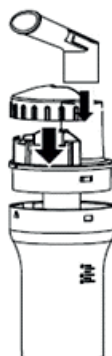SAFETY CONDITIONS

IMPORTANT INSTRUCTIONS ON SAFETY OF USE

PLEASE READ CAREFULLY AND KEEP FOR FUTURE REFERENCE

The warranty conditions are different, if the device is used for commercial purpose.

1. Before using the product please read carefully and always comply with the following instructions. The manufacturer is not responsible for any damages due to any misuse.
2. The product is only to be used indoors. Do not use the product for any purpose that is not compatible with its application.
3. The applicable voltage is 230V, ~50Hz. For safety reasons it is not appropriate to connect multiple devices to one power outlet.
4. Please be cautious when using around children. Do not let the children to play with the product. Do not let children or people who do not know the device to use it without supervision.
5. WARNING: This device may be used by children over 8 years of age and persons with reduced physical, sensory or mental capabilities, or persons without experience or knowledge of the device, only under the supervision of a person responsible for their safety, or if they were instructed on the safe use of the device and are aware of the dangers associated with its operation. Children should not play with the device. Cleaning and maintenance of the device should not be carried out by children, unless they are over 8 years of age and these activities are carried out under supervision.
6. After you are finished using the product always remember to gently remove the plug from the power outlet holding the outlet with your hand. Never pull the power cable!!!
7. Never leave the product connected to the power source without supervision. Even when use is interrupted for a short time, turn it off from the network, unplug the power.
8. Never put the power cable, the plug or the whole device into the water. Never expose the product to the atmospheric conditions such as direct sun light or rain, etc.. Never use the product in humid conditions.
9. Periodically check the power cable condition. If the power cable is damaged the product should be turned to a professional service location to be replaced in order to avoid hazardous situations.
10. Never use the product with a damaged power cable or if it was dropped or damaged in any other way or if it does not work properly. Do not try to repair the defected product yourself because it can lead to electric shock. Always turn the damaged device to a professional service location in order to repair it. All the repairs can be done only by authorized service professionals. The repair that was done incorrectly can cause hazardous situations for the user.
11. Never put the product on or close to the hot or warm surfaces or the kitchen appliances like the electric oven or gas burner.
12. Never use the product close to combustibles.
13. Do not let cord hang over edge of counter.
14. Before using the machine, please consult your physician to determine appropriate drug and dose.
15. Before first use and after each use, thoroughly wash the elements 2,3,4,6.
16. After washing and drying parts 2,3,4,6 store them in a clean place.
17. Do not shake or tilt the device when using.
18. Use only original accessories included in the set.

- 19 Do not use when only water was poured to the tank.
- 20 Do not wrap the power cord and air hose around the unit.
- 21 The product is not suitable for animals. For human use only.

DESCRIPTION (Fig.A)

1 On / Off switch	2 mask
3 Mouthpiece	4 air duct
5 Air Filter	6 nebulizer
7 Scale for measuring dose	8 Compressor

USE

- 1 Set the unit on a flat, stable surface.
- 2 Pour the metered medicament to the nebuliser reservoir (6).
- The appropriate dose can be measured by using the scale (7) (3ml min, max 8ml).
- 3 Put all the elements of the inhaler together (Fig. D) and make sure they are mounted properly.
- 4 Attach the mask (2) or the mouthpiece (3).
- NOTE: Make sure that the air hose (4) is properly attached (Fig.A)
- 5 Switch on the device by pressing button (1) and sit back while holding the nebulizer (6) with your hands. Place the mouthpiece / mask and breath deeply.
- 6 After the session, turn off and unplug the inhaler from the power source.
- 7 After the session, there is a possibility of condensation in the air duct - disconnect the air duct from the compressor (8) and nebulizer (6) and allow all elements to dry. Due to the air stream from the compressor, the air duct gets wet. This may cause mildew..

CLEANING

- 1 Before cleaning, disconnect the appliance from the power supply.
- 2 Dismount the elements.
- 3 Empty the nebulizer from the medicament residue.
- 4 Thoroughly wash parts 2,3,4,6 in hot water.
- 5 Use a wipe with a slightly damp cloth to clean the outer surfaces of the compressor.
- 6 Before reusing, thoroughly dry all parts of the device.
- 7 The device can be periodically wiped with a dry or slightly damp cloth.
- 8 Air filter should be changed every 60 days or when it changes color.

SPECIFICATIONS

Power supply: 230V ~ 50Hz

Power: 140W

Minimum / maximum capacity of the drug: 3ml / 8ml

Max pressure: 2.1-2.5bar

Pressure (nebulizer): 0.8-1.3bar

Compressor capacity: 6-8L / min

Minimum capacity: 0.2ml / min

% Of particles 0.1 - 10um

MMAD ~ 4µm

Noise level: <65dB (approx 1m)

Operating conditions: 10 ° C to 35 ° C, 15% to 85% RH

Storage: -20 ° C to + 65 ° C, 15% to 85% RH

Working / Storage 80KPa - 105KPa

Device is made in class II of insulation.

Device is compliant with EU directives:

- Low voltage directive (LVD)

- Electromagnetic compatibility (EMC)

Device marked CE mark on rating label

To protect your environment: please separate carton boxes and plastic bags and dispose them in corresponding waste bins. Used appliance should be delivered to the dedicated collecting points due to hazardous components, which may effect the environment. Do not dispose this appliance in the common waste bin.

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